

REMEDY® HIP SYSTEM

HIP MODULARITY

REMEDY® MODULAR STEMS


The REMEDY® Hip Spacer is part of the treatment foreseen in a two-stage procedure performed in the event of permanent prosthesis infection. The REMEDY® Hip Spacer implant is intended for temporary use only (180 days or less). It allows basic joint mobility and releases antibiotics into the joint area to protect the implant from bacterial colonization. A second surgery will be required at a later date to remove the REMEDY® Hip Spacer and replace it with a permanent hip joint implant.




Stainless steel rod reinforces the femoral stem

REMEDY® Hip Spacers:

- Single-use medical devices/ethylene oxide sterile
- Formed with bone cement (PMMA) and gentamicin

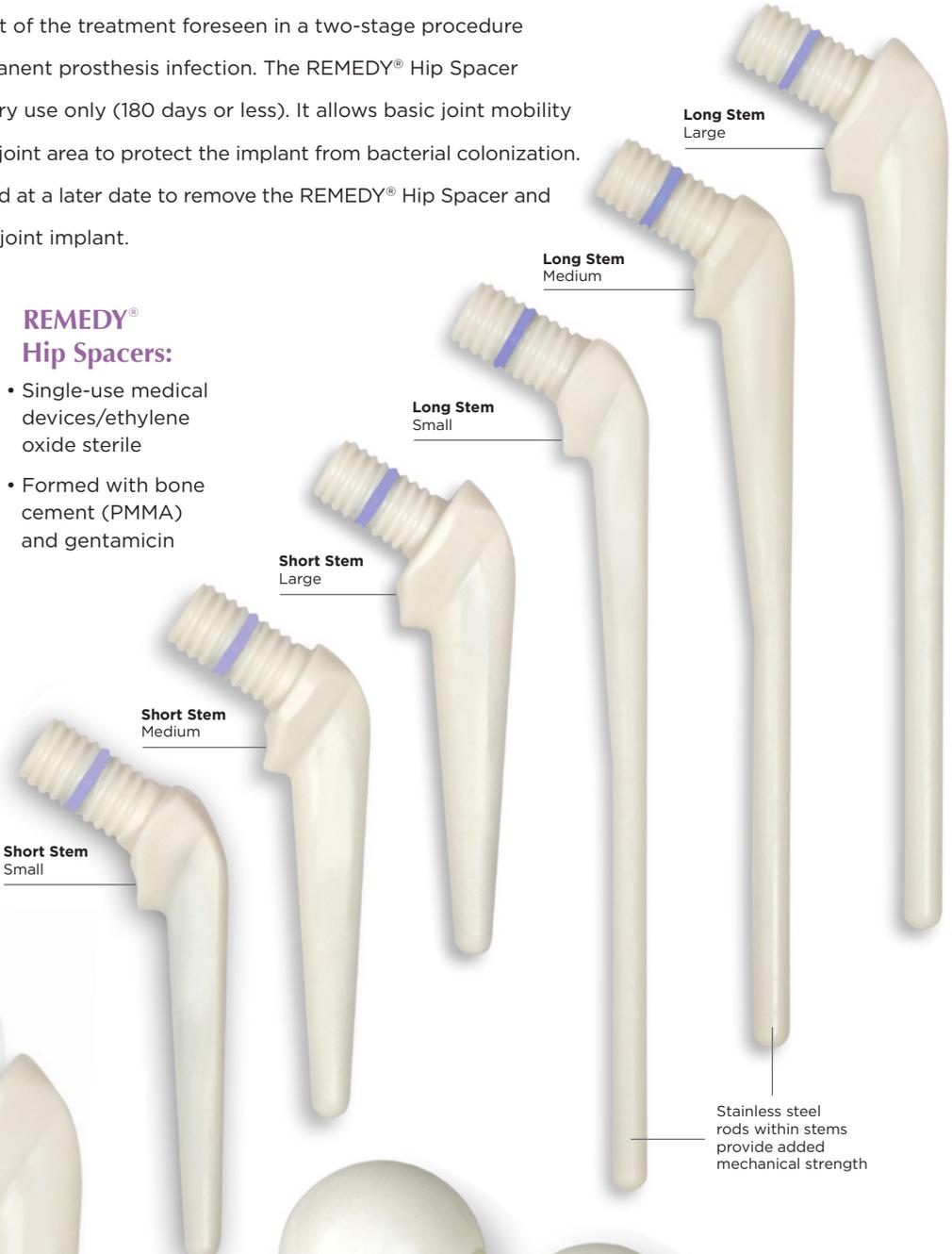

Stainless steel rods within stems provide added mechanical strength

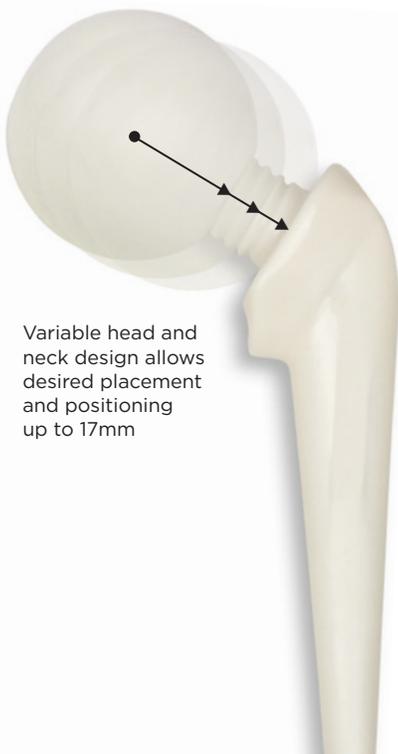

Variable head and neck design allows desired placement and positioning up to 17mm



60mm



54mm




46mm

REMEDY® MODULAR HEADS

Head sizes interchangeable with stems for surgical flexibility

HIP SPECIFICATIONS

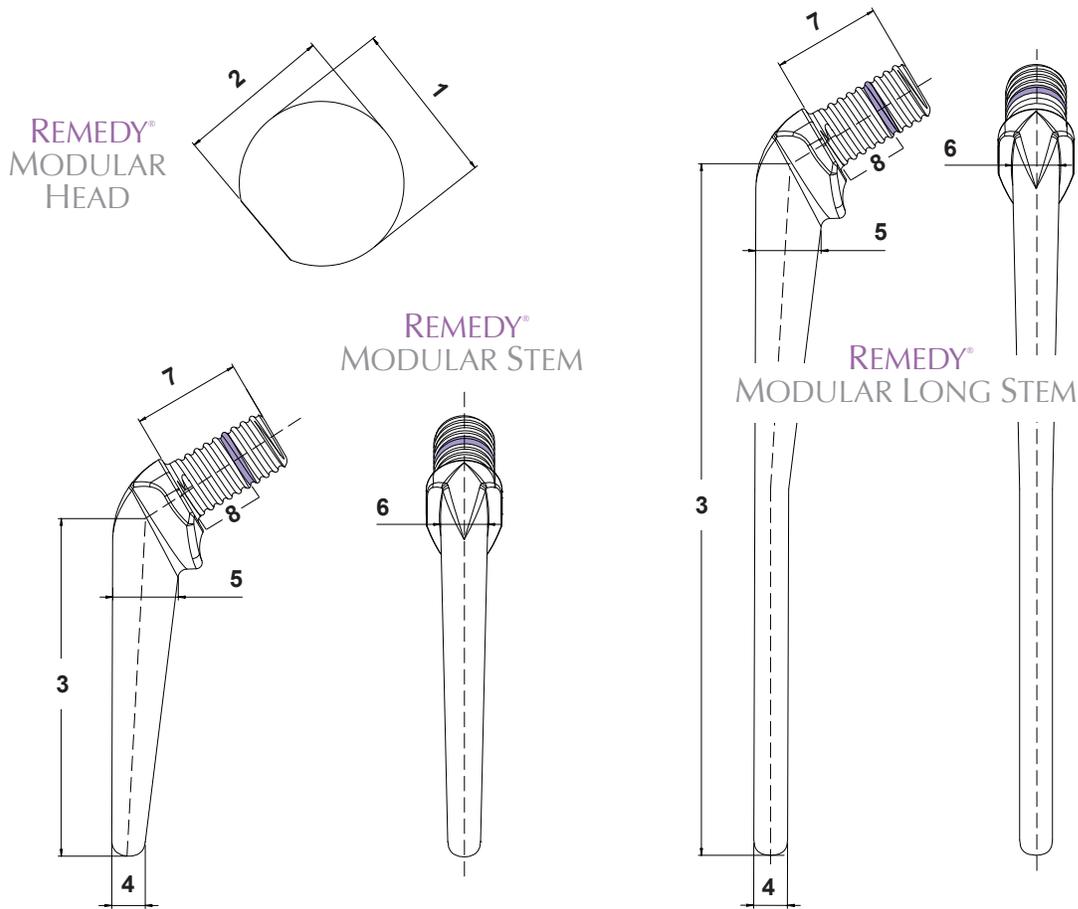

REMEDY® HIP SPACER

Description	Catalog #	(MM)								Gentamicin Base (g)
		1	2	3	4	5	6	7	8	
REMEDY® Modular Head 46mm	RHHDSM	46	42.3							0.9
REMEDY® Modular Head 54mm	RHHDMD	54	50.9							1.6
REMEDY® Modular Head 60mm	RHHDLG	60	57.3							2.3
REMEDY® Modular Stem - Small	RHSTSM			111	10	16.5	11.3	35.6	17	0.5
REMEDY® Modular Stem - Medium	RHSTMD			112	11	21.7	15.5	35.6	17	0.6
REMEDY® Modular Stem - Large	RHSTLG			117	11.5	24	16.5	35.6	17	0.7
REMEDY® Modular Long Stem - Small	RHLSSM			227	10	16.5	11.3	35.6	17	0.6
REMEDY® Modular Long Stem - Medium	RHLSMD			227	11	21.7	15.5	35.6	17	0.8
REMEDY® Modular Long Stem - Large	RHLSLG			231	11.5	24	16.5	35.6	17	0.9